Project #8 Microwave Based Fluid Bed Drying

Presenter: Doug Hausner

Faculty: Ben Glasser (Project Leader, Rutgers), Marianthi Ierapetritou (Rutgers)

IAB Meeting March 02-03, 2017











Problem Statement

- Current drying technologies being implemented for continuous solid dose manufacturing systems using conventional heated air drying do not operate on the same time scale (as other unit operations within the overall process.
- Drying using conventional methods is inherently slow due to its reliance on diffusion, which leads to characteristic response times that are not ideally suited for control of continuous manufacturing processes.
- Microwave drying is not new though, and has been utilized in industries outside of pharma for similar materials.
- ➤ In pharma, this technology may never have been fully explored for solid dose manufacturing applications because of the fear of product degradation. Today, while this concern may still remain, the process can be monitored and controlled much more tightly using PAT and feedback control.











Objectives

➤ Demonstrate the feasibility of utilizing microwave based drying modifications to traditional fluid bed drying unit operations that allows for drying times that more closely match the timescales of the other unit operations within a typical solid dose process.

Potential Future Directions if Funded for Year 2 and Beyond

Proof of concept implementation with integration of PAT and control









Materials/Methods

- A commercially available fluid bed dryer will be mounted with variable power magnetrons and aluminum insert
- Granulated compounds with well-known thermal properties and degradation products, including common pharmaceutical polymeric binders, will be chosen so that these aspects can be explored.
- DoEs will be performed on a range of selected materials to explore the effects of microwave intensity, duration, air-flow, and mass loading.
- A combination of Raman and Infra-red spectroscopy will be used offline to determine if any degradation products arise from any of the experimental points. Both types of spectroscopy will be used due to the fact that additional sensitivity will be available with Infra-red, but it would not be suitable for a PAT application.











Back of the Envelope

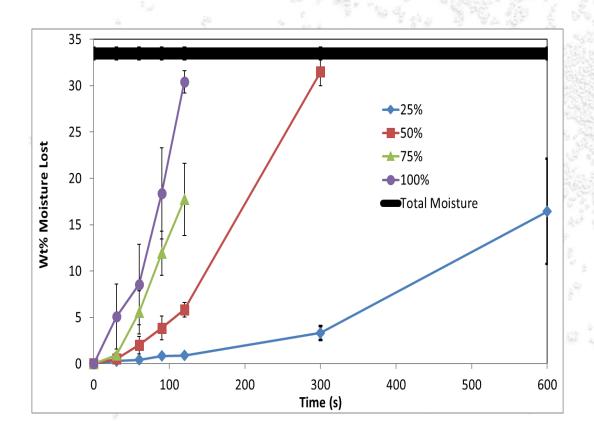


Figure of preliminary static microwave drying results of APAP-Lactose-PVC granules at 25-100% of 1250kW, no convection, suitable drying attained in 5 minutes or less











Anticipated Impact

➤ The development of continuous wet granulation processes with more efficient control strategies, ideally taking the response time from on the scale of 30 minutes to less than 5 minutes.









Q & A?









